


4 Research Court, Suite 300 Rockville, MD 20850 877-411-2041 Services@ibtbioservices.com

## Western Blot Detection

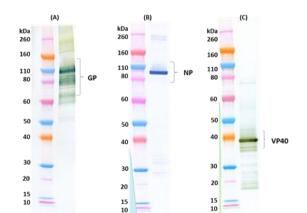

Western Blot was used to confirm the presence of EBOV-specific GP, NP and VP40 in the VLP. VLP sample was heat-denatured under reducing condition and loaded at  $1 \mu g/lane$ .

- (A) GP was detected as a broadly-diffused band 60-110 kDa using 50 ng/mL of anti-EBOV GP rabbit Pab (cat# 0301-015)
- (B) NP was detected as a single sharp band at  $\sim 90$  kDa using 100 ng/mL of anti-EBOV NP rabbit Pab (cat# 0301-012)
- (C) VP40 was detected as a single sharp band at  ${\sim}40~\rm kDa$  using 50 ng/mL of anti-EBOV VP40 rabbit Pab (cat# 0301-010)

## Ebola virus-like particles (EBOV VLP)

Catalog #: 0550-001

Lot #: 1810003

**Description:** Virus-like Particles (VLP's) expressing recombinant Ebola virus (EBOV) glycoprotein (GP), nucleoprotein (NP), and matrix protein (VP40). These VLP's are produced in Sf9 insect cells through infection with a recombinant baculovirus.

Storage: 2-3 weeks at -20°C, -80°C long term

**Size: 1 mg** of protein supplied in PBS (supplemented with arginine and glutamic acid) at a concentration of **2.403 mg/mL**.

Endotoxin: 2.30 Endotoxin Units/mg

**Relevance:** Since these VLP's mimic Ebola virus but do not contain genetic material, thus are not infectious and make an ideal candidate as a vaccine and as a tool to enhance filovirus research.

**Related Products:** IBT provides a wide array of anti-filovirus specific antibodies and other infectious disease reagents. Please see our website, www.ibtbioservices.com for more details.